



URBANSTREET

Orthosis Shoe. Different.






URBANSTREET

Style and function for everyday life.

The URBANSTREET orthosis shoe is not just new, but also different. It is the result of years of experience, which we gained during the construction of custom-made orthoses. The requirements and demands for an orthosis shoe of people wearing orthoses were the main focus.

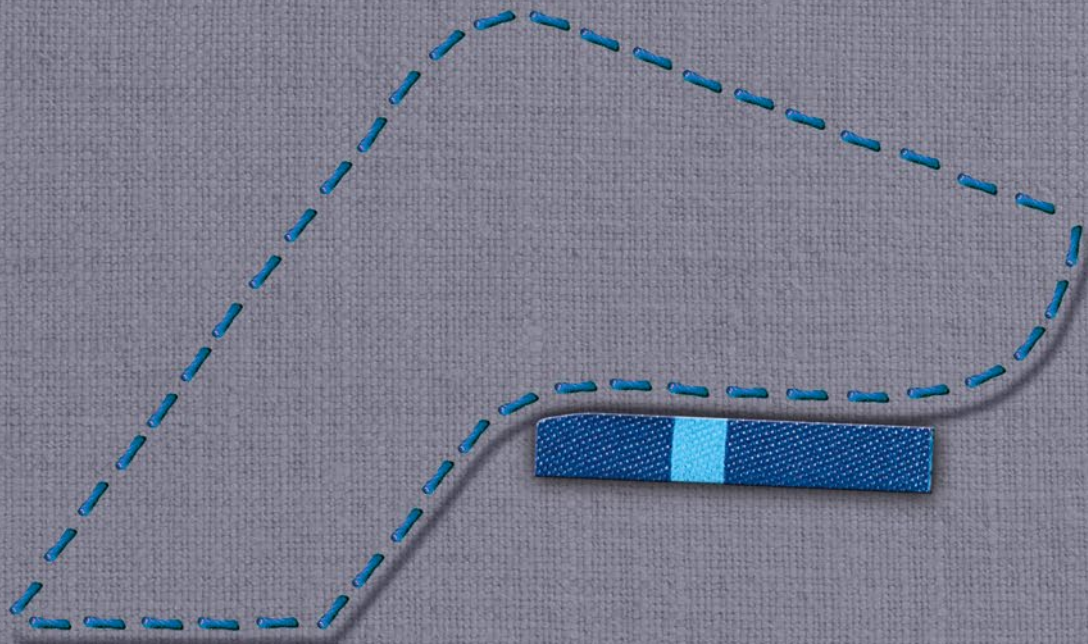


A unisex shoe that stands out from common orthosis shoes not only because of its extraordinary design features but as well because of its modern style. Casual through everyday life with URBANSTREET – the original with the boomerang.

URBANSTREET. Handmade with compelling details.



Heel Strap

The heel strap facilitates putting on and taking off the shoe.

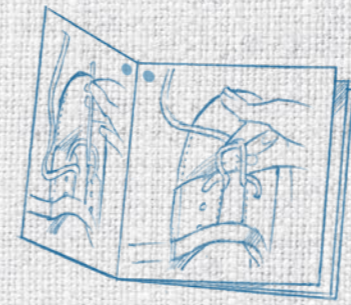

Guide for One-Handed Lacing

Each pair of shoes comes with a Guide for One-Handed Lacing as well as an additional lace for one-handed lacing, which offers support for patients with hemiplegia.



Individual Shortening and Wide Opening

Thanks to the elastic textile material and the combination of adjustable hook and loop fastener and laces, the shoe opening can be enlarged so that the shoe is comfortable to put on.

What Is Different?

The integrated elastic textile material offers comfortable space for orthoses with ankle joint and at the same time prevents the deformation of the shoe's upper.

Stability

The higher, anatomically formed upper and the integrated heel cap guarantee high stability.

Sturdy Fabric

The sturdy fabric made of densely woven cotton is durable, easy to clean and guarantees a carefree use in everyday life. If required, the shoes can be washed for 30 minutes in a washing machine at 30°C (86°F).



2 Insoles

The two insoles (5mm each) are removable and serve as volume compensation.

Multifunctional Outsole

The heelless rubber outsole insulates from coldness and wetness and has a slip-resistant, shock-absorbing tread.





Order Information

Art. No.	Insole Length
OS4921-32	209mm
OS4921-33	216mm
OS4921-34	223mm
OS4921-35	229mm
OS4921-36	235mm

Shoe Measurements		Use
pitch*	0mm	h-Cast
heel height**	14mm	Orthosis Configurator

Art. No.	Insole Length
OS4921-37	242mm
OS4921-38	249mm
OS4921-39	254mm
OS4921-40	261mm
OS4921-41	268mm
OS4921-42	275mm

Art. No.	Insole Length
OS4921-43	281mm
OS4921-44	287mm
OS4921-45	294mm
OS4921-46	301mm
OS4921-47	307mm
OS4921-48	314mm

* The height of the pitch corresponds to the height of the h-Cast plate to be used in the physical examination.

** The heel height is used in the Orthosis Configurator as part of the calculation basis.

only sold in pairs